

## Product datasheet

# Anti-Steroidogenic Factor 1/SF-1 antibody [EPR19744] ab217317

Recombinant RabMAb

[3 Images](#)

### Overview

<b>Product name</b>	Anti-Steroidogenic Factor 1/SF-1 antibody [EPR19744]
<b>Description</b>	Rabbit monoclonal [EPR19744] to Steroidogenic Factor 1/SF-1
<b>Host species</b>	Rabbit
<b>Tested applications</b>	<b>Suitable for:</b> IHC-P
<b>Species reactivity</b>	<b>Reacts with:</b> Human
<b>Immunogen</b>	Synthetic peptide within Human Steroidogenic Factor 1/SF-1 aa 200-300. The exact sequence is proprietary. Database link: <a href="#">Q13285</a>
<b>Positive control</b>	IHC-P: Human ovary, testis and adrenal cancer tissues.
<b>General notes</b>	Previously labelled as Steroidogenic Factor 1.  Our RabMAb <sup>®</sup> technology is a patented hybridoma-based technology for making rabbit monoclonal antibodies. For details on our patents, please refer to <a href="#">RabMAb<sup>®</sup> patents</a> .  This product is a <a href="#">recombinant rabbit monoclonal antibody</a> .

### Properties

<b>Form</b>	Liquid
<b>Storage instructions</b>	Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C long term. Avoid freeze / thaw cycle.
<b>Storage buffer</b>	Preservative: 0.01% Sodium azide Constituents: 59% PBS, 40% Glycerol, 0.05% BSA
<b>Purity</b>	Protein A purified
<b>Clonality</b>	Monoclonal
<b>Clone number</b>	EPR19744
<b>Isotype</b>	IgG

### Applications

## Applications

---

Our [Abpromise guarantee](#) covers the use of **ab217317** in the following tested applications.

The application notes include recommended starting dilutions; optimal dilutions/concentrations should be determined by the end user.

Application	Abreviews	Notes
IHC-P		1/2000. Perform heat mediated antigen retrieval with Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.

---

## Target

---

### Function

Transcriptional activator. Seems to be essential for sexual differentiation and formation of the primary steroidogenic tissues. Binds to the Ad4 site found in the promoter region of steroidogenic P450 genes such as CYP11A, CYP11B and CYP21B. Also regulates the AMH/Muellerian inhibiting substance gene as well as the AHCH and STAR genes. 5'-YCAAGGYC-3' and 5'-RRAGGTCA-3' are the consensus sequences for the recognition by NR5A1. The SFPQ-NONO-NR5A1 complex binds to the CYP17 promoter and regulates basal and cAMP-dependent transcriptional activity. Binds phosphatidylcholine (By similarity). Binds phospholipids with a phosphatidylinositol (PI) headgroup, in particular PI(3,4)P2 and PI(3,4,5)P3.

### Involvement in disease

Defects in NR5A1 are a cause of 46,XY disorder of sex development (46,XY DSD) [MIM:612965]; also known as XY sex reversal with or without adrenal failure. A congenital condition in which development of chromosomal, gonadal, or anatomic sex is atypical. 46,XY DSD is a disorder of gonadal (testicular) development, which may be complete or partial. The complete form includes streak gonads, normal mullerian structures, and normal female external genitalia. The partial form includes ambiguous external genitalia and partial development of mullerian and wolffian structures.

Defects in NR5A1 are a cause of adrenocortical insufficiency without ovarian defect (ACIWOD) [MIM:184757]. ACIWOD is characterized by severe 'slackness' muscular hypotonia. There is decreased sodium, increased potassium and elevated ACTH.

Defects in NR5A1 are the cause of premature ovarian failure type 7 (POF7) [MIM:612964]. An ovarian disorder defined as the cessation of ovarian function under the age of 40 years. It is characterized by oligomenorrhea or amenorrhea, in the presence of elevated levels of serum gonadotropins and low estradiol.

### Sequence similarities

Belongs to the nuclear hormone receptor family. NR5 subfamily.  
Contains 1 nuclear receptor DNA-binding domain.

### Post-translational modifications

Acetylation stimulates the transcriptional activity.


### Cellular localization

Nucleus.

---

## Images


---

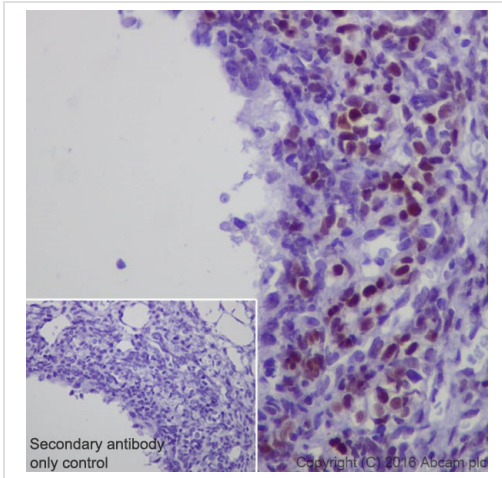

Immunohistochemical analysis of paraffin-embedded human ovary tissue labeling Steroidogenic Factor 1/SF-1 with ab217317 at 1/2000 dilution, followed by Goat Anti-Rabbit IgG H&L (HRP) Ready to use. Nuclear staining on human ovary is observed [PMID: 20660055] [PMID: 10370224]. Counter stained with Hematoxylin.

Secondary antibody only control: Used PBS instead of primary antibody, secondary antibody is Goat Anti-Rabbit IgG H&L (HRP) Ready to use.


Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-Steroidogenic Factor 1/SF-1 antibody [EPR19744] (ab217317)

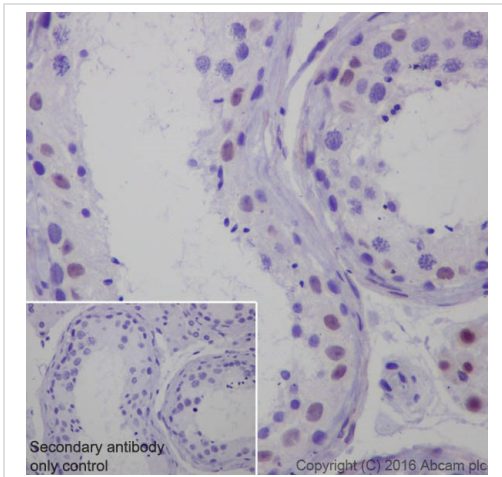

Immunohistochemical analysis of paraffin-embedded human testis tissue labeling Steroidogenic Factor 1/SF-1 with ab217317 at 1/2000 dilution, followed by Goat Anti-Rabbit IgG H&L (HRP) Ready to use. Nuclear staining on human testis is observed [PMID: 10370224]. Counter stained with Hematoxylin.

Secondary antibody only control: Used PBS instead of primary antibody, secondary antibody is Goat Anti-Rabbit IgG H&L (HRP) Ready to use.

Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-Steroidogenic Factor 1/SF-1 antibody [EPR19744] (ab217317)

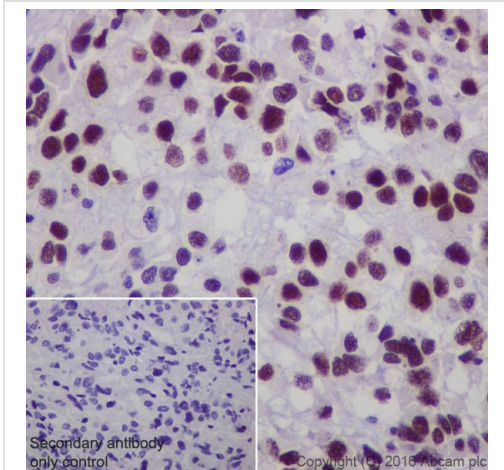

Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-Steroidogenic Factor 1/SF-1 antibody [EPR19744] (ab217317)

Immunohistochemical analysis of paraffin-embedded human adrenal cancer tissue labeling Steroidogenic Factor 1/SF-1 with ab217317 at 1/2000 dilution, followed by Goat Anti-Rabbit IgG H&L (HRP) Ready to use. Nuclear staining on human adrenal cancer is observed [PMID: 20660055]. Counter stained with Hematoxylin.

Secondary antibody only control: Used PBS instead of primary antibody, secondary antibody is Goat Anti-Rabbit IgG H&L (HRP) Ready to use.

**Please note:** All products are "FOR RESEARCH USE ONLY AND ARE NOT INTENDED FOR DIAGNOSTIC OR THERAPEUTIC USE"

### Our Abpromise to you: Quality guaranteed and expert technical support

---

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <https://www.abcam.com/abpromise> or contact our technical team.

### Terms and conditions

---

- Guarantee only valid for products bought direct from Abcam or one of our authorized distributors